## ENGLISH

### **Congratulations on purchasing your TriPollar® Stop X device!**

You have taken the first step toward flawless beautiful skin!  
This user manual describes the TriPollar® Stop X product and how to use it.

Supplementary information to be found at:  
<http://www.stop-age.com> (English website)  
<http://www.stop-age.cn> (Chinese website)

## TABLE OF CONTENTS

<b>INTRODUCING TRIPOLLAR® STOP X</b>	<b>8</b>
• OVERVIEW	
• THE TRIPOLLAR® STOP X KIT	
<b>USING TRIPOLLAR® STOP X</b>	<b>11</b>
• THE TRIPOLLAR® STOP X DEVICE	
• BEFORE YOU BEGIN	
• TREATING YOUR SKIN	
• RECOMMENDED TREATMENT GUIDELINES	
<b>TROUBLESHOOTING</b>	<b>20</b>
<b>CUSTOMER SUPPORT</b>	<b>20</b>
<b>TECHNICAL SPECIFICATIONS</b>	<b>21</b>
<b>APPENDIX</b>	<b>22</b>

## IMPORTANT NOTICE

© Pollogen Ltd. All rights reserved.

Design and technical specifications are subject to change without prior notice. Manufactured by Pollogen Ltd.

Patented technology.

The names Pollogen, TriPollar and Stop are registered trademarks of Pollogen Ltd.

All intellectual property rights in this publication are owned by Pollogen Ltd. and protected by applicable copyright laws and international treaty provisions. Pollogen Ltd. retains all rights not expressly granted. No part of this publication may be reproduced in any form whatsoever or used to make any derivative work without prior written approval by Pollogen Ltd.

All other trademarks are the property of their respective owners. Other company and brand products and service names are trademarks or registered trademarks of their respective holders.

## SAFETY

**Important:** Be sure to read the entire user manual prior to use of the product. Always adhere to all warnings and cautions. Be sure to follow all safety precautions below when using the TriPollar® Stop X product. To learn more about TriPollar technology, visit our English website at: [www.stop-age.com](http://www.stop-age.com) and Chinese website at: [www.stop-age.cn](http://www.stop-age.cn)

### ① WARNINGS AND CAUTIONS

TriPollar® Stop X is an electrical device that must be used properly. To reduce the risk of injury from electric shock, adhere to all warnings and cautions below:

- Do not touch the device should it fall into water. Unplug the device immediately with dry hands.
- Do not place or store the device near any source of water, such as a sink or bathtub.
- Do not submerge the device under water.
- When plugged in, do not leave the device unattended.
- Before use, ensure that the system and its power supply do not show any visible damage. Do not operate if the electric cord is damaged or if the system is disassembled.
- Only use the device with the electric cord and power supply provided.
- Keep the device out of the reach of children.
- Be sure that your hands are completely dry when operating the device.
- Turn the device OFF and then unplug the power supply from the electrical outlet after use or before cleaning and maintenance.

- Use TriPollar® Preparation Gel during treatment according to its instructions for use.

**SAFETY INSTRUCTIONS**

- Only use this device for its intended purpose, as described in this manual.
- Do not use the device over eyelids, closed eyes, thyroid area, lips, ears, breasts or genital areas.
- The sensation while using the device should be warm and pleasant. If you experience discomfort, discontinue use immediately.

Failure to follow the instructions in this manual may cause damage or burns to your skin.

**DO NOT USE THE DEVICE IF ANY OF THE FOLLOWING APPLY TO YOU:**

- Are under 18 years of age.
- Have a pacemaker or internal defibrillator.
- Have metal implants in the treatment area, such as gold threads. This does not apply to dental fillings or implants, which present no problem.
- Currently have cancer or a history of cancer, especially skin cancer or premalignant moles.
- Have severe concurrent conditions, such as cardiac disorders.
- Are pregnant or nursing.
- Have an impaired immune system due to immunosuppressive diseases such as AIDS, or are using immunosuppressive medications.
- Have any disease that is stimulated by heat, such as recurrent Herpes Simplex (within the treatment area).

- Have any active skin condition within the treatment area, such as sores, psoriasis, eczema or a rash.
- Have a history of keloid scarring, abnormal wound healing or fragile skin.
- Have had any surgical procedure, laser resurfacing or deep chemical peeling within the treatment area within the past three months or are still undergoing a process to complete healing.
- Have permanent fillers within the treatment area.  
If you have temporary dermal fillers within the treatment area, consult your physician before using the device.
- Have had Botox administered in the treatment area within the past three months.
- Have suspicious lesions within the treatment area.
- Currently use medications, herbal preparations, vitamins or food supplements that may cause fragile or dry skin.

## **POSSIBLE SIDE EFFECTS**

The following side effects are rare. However, should they occur, they are usually temporary. If any of the following conditions appear and last more than 24 hours, consult your physician:

- Pain
- Damage to natural skin texture (crusting, blistering, burns, fragile skin, bruising)
- Mild headache

For other aesthetic procedures or should you have any uncertainty regarding potential side effects, consult your physician before using the device.

## INTRODUCING TRIPOLLAR® STOP X

This chapter introduces the TriPollar® Stop X product and describes its components.

### OVERVIEW

#### **TriPollar® Stop X – The Next Generation of RF Home Treatments**

TriPollar® Stop X presents the latest in RF technology at the comfort of your home. Based on the groundbreaking TriPollar® Multi-RF technology, TriPollar® Stop X affects all layers of the dermis homogeneously for faster, safer and more effective treatment. It delivers professional-level aesthetic results that are both immediately visible and long-lasting. TriPollar® Stop X users enjoy the combined benefits of a unique non-conductive gel that ensures the delivery of focused RF energy into the right layers of the dermis using the newly designed electrodes. A proprietary and innovative temperature sensor ensures an effective and safe treatment while stimulating dermal activity that accelerates the natural production of collagen, the protein responsible for the skin's smooth and supple texture. Using TriPollar® Stop X regularly according to its treatment protocol can effectively postpone the signs of aging.

## THE TRIPOLLAR® STOP X KIT


The TriPollar® Stop X kit contains the following components:

1. TRIPOLLAR® STOP X KIT (FIG.1)

FIG.1

2. TRIPOLLAR® STOP X DEVICE (FIG.2)
3. TRIPOLLAR® STOP X POWER SUPPLY AND ELECTRIC CORD (FIG.3)
4. TRIPOLLAR® STOP X ELECTRICAL OUTLET ADAPTERS (FIG.4)
5. TRIPOLLAR® PREPARATION GEL (FIG.5)  
After the container has been opened, the gel can be used for up to 12 months.  
**Important:** It is recommended to perform a patch test 24 hours prior to initial use. See page 13 for details.
6. TRIPOLLAR® STOP X USER MANUAL (FIG.6)
7. TRIPOLLAR® STOP X WARRANTY POLICY (FIG.6)
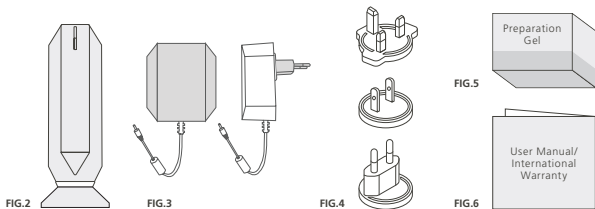
## USING TRIPOLLAR® STOP X

This chapter describes the TriPollar® Stop X device and how to use it.

### THE TRIPOLLAR® STOP X DEVICE

TriPollar® Stop X is a facial renewal and rejuvenation device. The device is used to increase and accelerate collagen production in your skin for rejuvenation effect.

TriPollar® Stop X focuses heat deep within the dermis of your skin, as shown below.


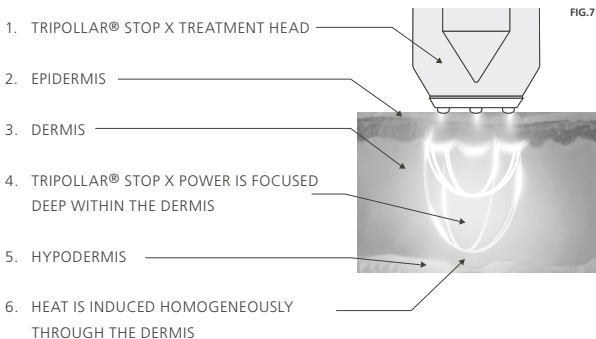
FIG.8 Shows the components comprising the TriPollar® Stop X device.

1. DEVICE BODY
2. TREATMENT HEAD
3. ELECTRODE
4. TEMPERATURE SENSOR
5. CONTROL PANEL
6. POWER SUPPLY CONNECTION
7. GREEN INDICATOR LIGHT
8. ORANGE INDICATOR LIGHT
9. ON/OFF BUTTON
10. LOW RF
11. HIGH RF
12. PROTECTIVE COVER
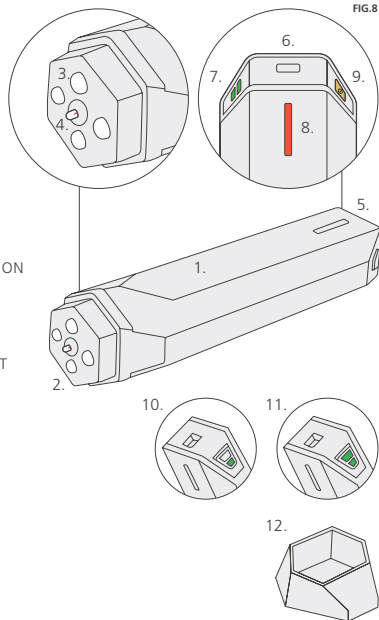
## BEFORE YOU BEGIN

Before using the TriPollar® Stop X device, do the following.

1. Open the box containing the device, and remove the device and its power supply.
2. Insert the adapter onto the electrical posts (Fig.9)
3. Push the adapter into the slot of the PS until it securely "clicks into place" and cannot be removed easily (Fig.10)
4. Connect the power supply to the device (Fig.11)
5. Plug the power supply into an electrical outlet (Fig.12)
6. Remove the TriPollar® Preparation Gel from its box.  
**Important:** Before using the device for the first time, it is recommended to perform a patch test 24 hours prior to initial use. To do so, apply a small amount of TriPollar® Preparation Gel to the inner arm and treat the area with the TriPollar® Stop X device. Ensure that the skin does not have a negative reaction to the gel or treatment.
7. Remove the protective cover from the device.

**You are now ready to treat your skin using TriPollar® Stop X, as described on page 15.**

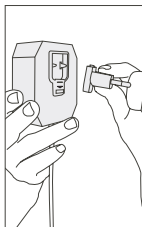
FIG.9

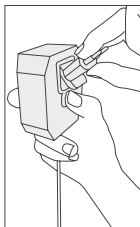
FIG.10

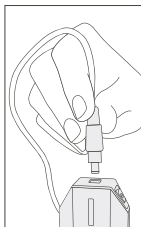
FIG.11

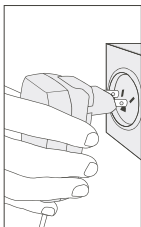
FIG.12

## TREATING YOUR SKIN

To treat your skin using the TriPollar® Stop X device:

1. Cleanse and dry the treatment area.
2. Apply a very thin, even layer of TriPollar® Preparation Gel over the entire treatment area (face and/or neck). This gel must be used during the entire treatment session. Your skin may feel a warm sensation when the gel is applied.
3. Turn the device ON by pressing the Power button (Fig.13). Make sure that the green indicator light turns ON.  
**Tip:** A long press of the Power button turns the device OFF.
4. Select the power setting (Low/High) that is most comfortable for you. When the green indicator light shows a single green light illuminated (Fig.14), the device is operating at the Low power level. You can press the Power button again to change to the High power level (Fig.15).

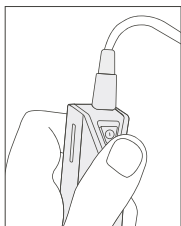
FIG.13

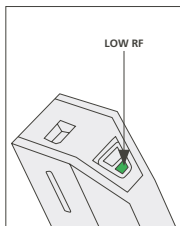
FIG.14

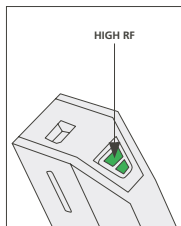
FIG.15

**Important:** During your first treatment session, we recommend that you use the Low power setting and then adjust the power level according to your skin's tolerance for subsequent treatments.

Treatment:

5. Divide the face to 3 optional treatment areas (Fig.16):
  - Right cheek
  - Left cheek
  - Forehead
6. Place the device on your skin, make sure to apply a slight pressure in order to assure a full contact between the electrodes and the skin.
7. Start moving the device in a very slow circular motion on a small part of your treatment area (approximately twice the size of the electrodes), slowly heating your skin to the optimal temperature (Fig.17).

**Important:** Be sure to constantly move the device slowly over the skin to ensure proper heating of the treated area. Do not leave the device on one place without movement.

8. Once the orange indicator light illuminates followed by a 'beep' sound, you have reached the optimal temperature. Continue the same slow circular motion for a few moments and then gradually enlarge your movements in the treatment area.  
Proceed with this technique covering the entire treatment area, making sure the orange indicator light is illuminated as much as possible during the treatment.  
As you move from one treatment area to another (Fig.16), the orange indicator light may turn-off if the temperature of the skin is not warm enough.
9. Continue working on the entire treatment area in a slow circular motion **for 4-5 minutes**. Make sure to maintain a uniform heating while covering the entire treatment area. Treatment can be augmented by "lifting" movements from the jaw line upwards toward your cheeks (Fig.18).
10. Once you finished the 4-5 minutes of treatment, you can proceed to the next treatment area repeating steps 6-9.

FIG. 16

FIG. 17

FIG. 18

11. As a natural result of the treatment, the skin may show some redness (erythema) and minor skin swelling (edema) in the treatment area. Both should disappear after a short time.
12. Press the Power button for two seconds to turn the device OFF.
13. Unplug the power supply from the electrical outlet.
14. Gently remove the TriPollar® Preparation Gel from your skin and wash your face. You can apply a moisturizer if desired.
15. Clean the electrodes of the device with a damp cloth or paper towel and cover the device tip with its protective cover.

**Warning: Do not submerge the TriPollar® Stop X device under water.**

**Safety Mechanism:**

The device's automatic temperature-optimization mechanism protects the skin from overheating. When the skin has reached the maximal temperature, the energy source will shut down automatically. The indication is - green LED off.

The device automatically turns the energy source back ON, when the skin temperature drops below the maximal temperature or when you move the device to an adjacent treatment area, where the skin temperature is lower. The indication is - green LED on.

**Important:**

The orange light is ON when the temperature of the skin in the treatment area is sufficiently heated and is OFF when it is not.

## RECOMMENDED TREATMENT GUIDELINES

The following treatment guidelines are recommended for using the TriPollar® Stop X device:

- Treating the entire face takes 12-15 minutes. For shorter treatment times focus on selected areas of concern. Treatment times differ for each person.
- For optimal results, use the TriPollar® Stop X device 2-3 times a week for a period of 6-8 weeks.
- Make sure that you allow at least one day in between treatments on the same treatment area.
- To maintain long-term results after finishing the full course of treatments, we recommend that you use your TriPollar® Stop X device 2-4 times a month or as needed.
- If the orange indicator light and sound indicators do not appear within five minutes, or if you do not feel any heating, try moving the device as slowly as possible without stopping the motion. Slow, circular, constant movements are more effective than rapid movements. The slower the movement, the faster the treatment area heats internally.
- If the orange indicator light and sound indicators turn ON immediately or if you feel excessive heat, select a lower power level.
- If the orange indicator light and sound indicators do not appear after 3–5 minutes, select a higher power level. If the indicators still do not appear after increasing to the High power setting, try moving the device more slowly, without stopping the motion.

## TROUBLESHOOTING

This chapter describes how to troubleshoot problems that may arise when using the TriPollar® Stop X device.

**If your TriPollar® Stop X device is not working, take the following steps before contacting the service center:**

- Ensure that you have removed the device's protective cover.
- Make sure that the electric cord is securely connected to the device and the Power supply adapter is firmly connected.
- Press the Power button and verify that the green indicator light illuminates.
- Ensure that you applied TriPollar® Preparation Gel according to instructions described in this manual. The Tripollar® Stop X device does not work properly when the preparation gel is not applied.

## CUSTOMER SUPPORT

This chapter provides customer support contact information for the Tripollar® Stop X device.

For service, refer to the store where you purchased the TriPollar® Stop X device or visit our website at: [www.stop-age.com](http://www.stop-age.com) and Chinese website at: [www.stop-age.cn](http://www.stop-age.cn) When contacting a local service center, have the device serial number ready. The serial number is located on the top panel of the device. Customer satisfaction is our top priority. To help us provide you with the best support possible, please send us your comments and suggestions. Contact information is provided on the back cover of this manual, or you can contact us via our website (see above).

## TECHNICAL SPECIFICATIONS

This chapter presents the specifications for the TriPollar® Stop X system.

TABLE 1

### Technical Specifications

Weight:	85gm
Physical Dimensions:	H=134mm L=51mm W=32mm
TriPollar RF Frequency:	1.0MHz-1.25MHz

TABLE 2

### Operating conditions

Temperature:	+5°C to +30°C
Humidity:	20RH to 80RH
Rating:	100~240V AC, 0.6A, 50/60Hz DC 8V 1.5A

TABLE 3

### Transport and Storage

Temperature:	0°C to 40°C
Humidity:	10RH to 90RH

TABLE 4

### Electrical Safety

Power supply is classified as Class II double insulation

## APPENDIX

This appendix provides the hazardous materials declaration for the TriPollar® Stop X device.

### Restriction of Hazardous Substance Material Declaration

TABLE 5 - IS IN ACCORDANCE WITH THE PROVISIONS OF SJ/T 11364

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- x Indicates that the named hazardous substance contained in at least one of the homogeneous materials used for this component may be above the limit requirement of GB/T 26572. The restricted materials are used in accordance with manufacturing needs. Use of such materials has no impact on the device's safety and no harmful effect on device users. **This declaration is included to remind device users that the device must be disposed of reasonably and in accordance with the environmental policy of your locale.**
- Indicates that the named hazardous substance contained in all of the homogeneous materials for this component is below the limit requirement of GB/T 26572.

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